

NEW!



Grooved Mushroom Joystick Handle

A revolutionary new joystick handle, the Grooved Mushroom takes the best features of several previous Bodypoint® models and advances them to the next level. Its hand-filling shape, top and side grooves, and height-adapability make it ultra versatile to meet a wide range of user needs.



Top and side grooves give superior grip and control with limited hand function.

Textured, rounded top is comfortable and secure for users grasping from above.

Perfect height for upright gripping, yet easily trimmed shorter for palmar hand position or fingertip control. (See cutting instructions below.)

Grooved shaft and ability to tuck fingers below the knob make driving easy for users grasping from the side.

Springy and tough polyurethane is safe (non-cytotoxic), easy to grip, and a breeze to clean by hand or in dishwasher.


Flexible shaft dampens vibration in driving and reduces risk of damage to electronic controls.

Part Numbers

PC110A For 4.8 mm (3/16") stem size

PC110B For 6.4 mm (1/4") stem size

Easily cut to desired height with a utility knife

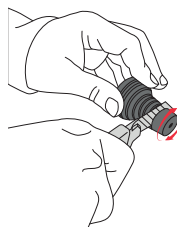

60 mm
(2-3/8")

28 mm
(1-1/8")